




Pressure Redistribution Mattress

NV-PPM9



Heel Slope

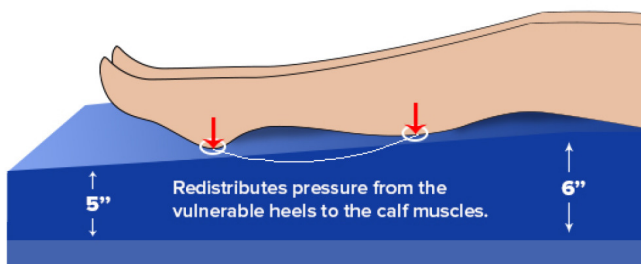

Model Numbers

NV-PPM9-80	36"W x 80"L x 6"H	Without Raised Side
NV-PPM9-80-RS	36"W x 80"L x 6"H	With Raised Side
NV-PPM9-84	36"W x 84"L x 6"H	Without Raised Side
NV-PPM9-80-42	42"W x 80"L x 6"H	Without Raised Side
NV-PPM9-84-42	42"W x 84"L x 6"H	Without Raised Side

Features:

- High quality therapeutic foam mattress provides pressure redistribution and shear/friction reduction
- Premium therapeutic foam mattress maximizes value and comfort with multiple layers of latex free foam
- Cover material is Stretch, Fluid-Resistant Nylon Top | Waterproof Vinyl Bottom
- Full top layer of open cell, visco-elastic memory foam improves immersion and pressure redistribution
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- High density foam perimeters provide solid edge support and facilitate safe patient transfers
- Available with or without raised side rails
- Heel section is sloped with channel cut foam center offloading pressure from vulnerable heels
- Removable multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Meets flammability standard CFR 1632 and CFR 1633
- 450 lbs weight capacity